**Reaktive Inhaltsstoffe:** Ein Test enthält pro cm<sup>2</sup>: **Leukozyten:** Indoxylester 15,5 µg, Methoxy-morpholinobenzodiazonium Salz 5,5 µg; **Nitrit:** Hydroxytetrahydrobenzo-quinolin 33,5 µg, Sulfanilamid 29,1 µg; **pH:** Bromthymolblau 13,9 µg, Methylrot 1,2 µg, Phenolphthalein 8,6 µg; **Protein:** Tetrachlorphenoltetrabromosulfophthalein 13,9 µg; **Glucose:** Tetramethylbenzidin 103,5 µg, GOD 6 U, POD 35 U; **Keton:** Nitroprussid-natrium 1572 µg, Glycin 4,2 mg; **Urobilinogen:** Methoxybenzodiazonium Salz 67,7 µg; **Bilirubin:** Dichlorbenzodiazonium Salz 16,7 µg; **Blut:** Tetramethylbenzidin 52,8 µg, Dimethylidihydroperoxyhexan 2972 µg.

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CE

COMBUR-TEST and COBAS are trademarks of Roche.

Roche Diagnostics GmbH  
D-68298 Mannheim, Germany  
www.diavant.com

Roche

0 4907957001(01)  
R1 (black)

Nine-patch test strip for the semiquantitative determination of leukocytes, nitrite, pH, protein, glucose, ketone bodies, urobilinogen, bilirubin, and blood in urine by visual reading

For professional use

FOR IN VITRO DIAGNOSTIC USE

Additionally required materials: A vessel for collection of urine

Instructions for use:

- Use fresh urine that has not been centrifuged. Thoroughly mix the urine sample. The sample should be at room temperature when the test is performed and should not have been standing for more than 2 hours.
- Take a test strip out of the container. Close the container again with the original desiccant stopper immediately after removal of the strip. This is important as otherwise the test areas may become discolored due to moisture and incorrect results may be obtained.
- Briefly (about 1 second) dip the test strip into the urine making sure that all test areas are moistened.
- When withdrawing the test strip, wipe the edge against the rim of the vessel to remove excess urine.
- After 60 seconds (60–120 seconds for the leukocyte test area), compare the reaction colors of the test areas with the colors on the label and assign always the value of the nearest color block. Compare the 9th (blood) test area with both color scales as separate color scales are given for erythrocytes and hemoglobin. Any color charges appearing only along the edges of the test areas, or developing after more than 2 minutes, do not have any diagnostic significance.

Test principles and notes on individual parameters

**Leukocytes:** The reaction detects the presence of esterases that occur in granulocytes. These enzymes cleave an indoxyl ester, and the indoxyl so liberated reacts with a diazonium salt to produce a violet dye. Both **intact** and **lysed** leukocytes are detected.

Reaction colors that cannot be unequivocally classed as negative or approx. 10–25 WBCs/ $\mu$ L urine after 60 seconds can generally be more easily assessed after 120 seconds.

The reaction is not affected by bacteria, trichomonads or erythrocytes present in the urine. Formaldehyde (stabilizer) and medication with antibiotics containing imipenem, meropenem or clavulanic acid may cause false-positive reactions. If the urine specimen is strongly colored (for example due to the presence of bilirubin or nitrofurantoin), the reaction color may be masked. Urinary protein excretion in excess of 500 mg/dL and urinary glucose excretion in excess of 1 g/dL may diminish the intensity of the reaction color, as can medication with antibiotics containing cephalixin or gentamicin if administered in high daily doses.

**Nitrite:** The most common organisms causing urinary tract infections, E.coli and most urinary pathogens, convert dietary nitrate to nitrite, which produces a pink to red coloration of the test area. The reaction thus indirectly detects the presence of nitrite-forming organisms in the urine. Even a slight pink coloration is an indication of significant bacteriuria.

Prolonged urinary retention in the bladder (4–8 hrs; ideally over night) is essential for a valid result. Administration of antibiotics or chemotherapeutics should be discontinued 3 days before the test. The test is based on the principle of Griess' test and is specific for nitrite.

**pH:** The pH values found most frequently in fresh urine from healthy individuals lie between 5 and 6. The test is specific for the detection of hydronium ions, the pH being the negative common logarithm of the hydronium ion concentration. The test pad contains the indicators methyl red, phenolphthalein and bromothymol blue.

**Protein:** The test is based on the principle of the protein error of pH indicators and is particularly sensitive to albumin.

Quinine, quinidine, chloroquine and tolbutamide do not affect the test, nor does a high pH (up to pH 9). False-positive results may be obtained after infusion of polyvinylpyrrolidone (blood substitute), or if the urine specimen collection vessel contains residues of disinfectants based on quaternary ammonium compounds or chlorhexidine.

**Glucose:** The determination of glucose is based on the specific glucose-oxidase/ peroxidase reaction.

This test is independent of the pH and specific gravity of the urine and is not affected by the presence of ketone bodies. The effect of ascorbic acid (vitamin C) has been largely eliminated so that at glucose concentrations of 100 mg/dL (5.5 mmol/L) and above even high ascorbic acid concentrations are not likely to give false-negative results.


**Ketone bodies:** This test is based on the principle of Legal's test and is more sensitive to acetoacetic acid than to acetone. Phenylketones and phthalein compounds produce red colors on the test area. These are, however, quite different from the violet colors produced by ketone bodies. Captopril, mesna (2-mercaptoethanesulfonic acid sodium salt) and other substances containing sulphydryl groups may produce false-positive results.

**Urobilinogen:** A stable diazonium salt reacts almost immediately with urobilinogen to give a red azo dye. No discoloration of the test area or colors lighter than that shown for 1 mg/dL (17  $\mu$ mol/L) constitute a normal finding. The test is specific for urobilinogen and is not susceptible to the interfering factors known to affect Ehrlich's test. Larger amounts of bilirubin produce a momentary yellow coloration of the test area which may turn green to blue after about 60 seconds.

**Bilirubin:** The test for bilirubin is based on the coupling of bilirubin with a diazonium salt to give an azo dye. Even the slightest pink coloration constitutes a positive, i.e. pathologic, result. Other urinary constituents produce a more or less intense yellow discoloration.

**Blood:** Hemoglobin and myoglobin catalyze the oxidation of the indicator by an organic hydroperoxide contained in the test paper. Separate color scales for erythrocytes and hemoglobin are given on the label of the test strip container. Individual to closely packed green dots on the yellow test area are indicative of intact erythrocytes. Hemoglobin, hemolyzed erythrocytes, and myoglobin are indicated by a uniform green coloration of the test area.

Ascorbic acid (vitamin C) has virtually no influence on test results.

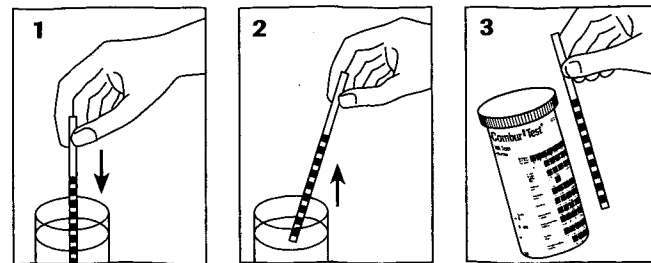**Note:** In women the test for blood may be falsified from 3 days before to 3 days after a period. It is therefore advisable not to perform the test during this time. After physical activity, e.g. strenuous jogging, raised values for erythrocytes and protein may occur without being signs of disease.

| Parameter            | Visual Reading                                   |                                                                                                            |                                                                                                  |
|----------------------|--------------------------------------------------|------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------|
|                      | Range                                            | Practical detection limit                                                                                  | Accuracy                                                                                         |
| Leukocytes           | Negative - approx. 500 WBC/ $\mu$ L (3+)         | 10–25 WBC/ $\mu$ L                                                                                         | $\geq$ 90% compared with counting chamber                                                        |
| Nitrite              | Negative - positive (1+)                         | 0.05 mg/dL (11 $\mu$ mol/L)                                                                                | $\geq$ 90% for 10 <sup>7</sup> gram-positive organisms compared with Griess' test                |
| pH                   | 5–9                                              |                                                                                                            | $\geq$ 95% compared with pH-meter                                                                |
| Protein              | Negative - 500 mg/dL (5 g/L; 3+)                 | 6 mg albumin/dL                                                                                            | 90% compared with radial immunodiffusion                                                         |
| Glucose              | Normal - 1000 mg/dL (55 mmol/L; 4+)              | 40 mg/dL (2.2 mmol/L)                                                                                      | $\geq$ 90% compared with hexokinase method                                                       |
| Ketone bodies        | Negative - 150 mg/dL (15 mmol/L; 3+)             | For acetoacetic acid 5 mg/dL (0.5 mmol/L)                                                                  | $\geq$ 85% compared with photometric enzymatic determination of acetate                          |
| Urobilinogen         | Normal - 12 mg/dL (200 $\mu$ mol/L; 4+)          | 0.4 mg/dL (7 $\mu$ mol/L)                                                                                  | $\geq$ 95% compared with Watson & Henry method                                                   |
| Bilirubin            | Negative - approx. 6 mg/dL (100 $\mu$ mol/L; 3+) | 0.5 mg/dL (9 $\mu$ mol/L)                                                                                  | $\geq$ 85% compared with total bilirubin determination by Jendrassik's method (direct bilirubin) |
| Blood and hemoglobin | Negative - approx. 250 RBC/ $\mu$ L (4+)         | Intact erythrocytes: 5 RBC/ $\mu$ L Hemoglobin or hemolysed erythrocytes: corresponding to 10 RBC/ $\mu$ L | $\geq$ 90% compared with counting chamber                                                        |

**Reactive components per cm<sup>2</sup>:** **Leukocytes:** indoxyl ester 15.5  $\mu$ g, methoxy-morpholinobenzene diazonium salt 5.5  $\mu$ g; **Nitrite:** hydroxytetrahydrobenzoquinoline 33.5  $\mu$ g, sulfanilamide 29.1  $\mu$ g; **pH:** bromothymol blue 13.9  $\mu$ g, methyl red 1.2  $\mu$ g, phenolphthalein 8.6  $\mu$ g; **Protein:** tetrachlorophenoltetrabromosulfophthalein 13.9  $\mu$ g; **Glucose:** tetramethylbenzidine 103.5  $\mu$ g, GOD 6 U, POD 35 U; **Ketone bodies:** nitroprusside sodium 1572  $\mu$ g, glycine 4.2 mg; **Urobilinogen:** methoxybenzene diazonium salt 67.7  $\mu$ g; **Bilirubin:** dichlorobenzene diazonium salt 16.7  $\mu$ g; **Blood:** tetramethylbenzidine 52.8  $\mu$ g, dimethyldihydroperoxyhexane 297.2  $\mu$ g.

**Please note:** Diagnosis or therapy should never be based on one test result alone but should be established in the context of all other medical findings. Not all effects of drugs or their metabolites upon the individual tests are known. In doubtful cases, it is therefore advisable to repeat the test after discontinuation of the medication. Use only clean, well-rinsed vessels to collect urine. False-positive readings, particularly for glucose, protein and blood, can result from residues of detergent or strongly oxidizing disinfectants in the specimen collection vessel.

Do not add stabilizers to the urine. Do not expose urine specimens to sunlight as this induces oxidation of bilirubin and urobilinogen and hence leads to artificially low results for these two parameters. Drugs that turn red in an acid environment (e.g. phenazopyridine) may produce false-positive readings or reddish discolorations on the test areas for nitrite, protein, urobilinogen, and bilirubin. Large quantities of ascorbic acid (vitamin C) can lead to low or false-negative results for nitrite and bilirubin.

**Storage and shelf-life:** Do not store the Combur<sup>9</sup>Test Pack at temperatures below +2°C or above +30°C. In the original container the test strips are stable up to the date printed on the pack, even the container has been opened.

**Disposal:** Please dispose of used test strips according to the safety regulations applicable at your facility.

The stopper of the test strip container is filled with a non-toxic silicate-based desiccant. If inadvertently ingested it should be flushed down with plenty of water.

**Presentation:** Packs of 50 test strips (REF 04510038191), 100 (REF 04510046040).

For an explanation of the symbols used and a list of references please refer to the end of this insert.

Last updated: 2006-10

Distributed by:

Roche Diagnostics Ltd  
Charles Avenue, Burgess Hill, RH15 9RY, United Kingdom

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